

Approach to mobile EEG sensor validation through resting state simultaneous data collection

Ekaterina Kutafina^{Corresp., 1, 2}, **Yannic Titgemeyer**¹, **Alexander Brenner**¹, **Rainer Surges**³, **Stephan M Jonas**⁴

¹ Department of Medical Informatics, Uniklinik RWTH Aachen, Aachen, Germany

² Faculty of Applied Mathematics, AGH University of Science and Technology, Krakow, Poland

³ Department of Epileptology, University Hospital of Bonn, Bonn, Germany

⁴ Department of Informatics, Technical University of Munich, Garching, Germany

Corresponding Author: Ekaterina Kutafina
Email address: ekutafina@mi.rwth-aachen.de

Development of mobile sensors brings new opportunities to medical research. In particular, mobile electroencephalography (EEG) devices can be potentially used in low cost screening for epilepsy and other neurological and psychiatric disorders. The necessary condition for such applications is thoughtful validation in the specific medical context. As part of validation and quality assurance, we develop a computer-based analysis pipeline, which aims to compare the EEG signal acquired by a mobile EEG device to the one collected by a medically approved clinical-grade EEG device. Both signals are recorded simultaneously during 30 minutes long sessions in resting state. The data are collected from 22 patients with epileptiform abnormalities in EEG. In order to compare two multichannel EEG signals with differently placed references and electrodes, a novel data processing pipeline is proposed. It allows deriving matching pairs of time series which are suitable for similarity assessment through Pearson correlation. The average correlation of 0.55 is achieved on a test dataset, which can be considered a promising result, taking the positions shift due to the simultaneous electrode placement into account.

Approach to mobile EEG sensor validation through resting state simultaneous data collection.

Ekaterina Kutafina^{1,2}, Yannic Titgemeyer¹, Alexander Brenner¹, Rainer Surges³, Stephan Michael Jonas⁴

¹Department of Medical Informatics, Uniklinik RWTH Aachen, Aachen, Germany

² Faculty of Applied Mathematics, AGH University of Science and Technology, Krakow, Poland

³Department of Epileptology, University Hospital of Bonn, Bonn, Germany

⁴Department of Informatics, Technical University of Munich, Garching, Germany

Corresponding Author:

Ekaterina Kutafina

Pauwelsstrasse 30, 52057 Aachen, Germany

email: ekutafina@mi.rwth-aachen.de

Abstract: Development of mobile sensors brings new opportunities to medical research. In particular, mobile electroencephalography (EEG) devices can be potentially used in low cost screening for epilepsy and other neurological and psychiatric disorders. The necessary condition for such applications is thoughtful validation in the specific medical context. As part of validation and quality assurance, we develop a computer-based analysis pipeline, which aims to compare the EEG signal acquired by a mobile EEG device to the one collected by a medically approved clinical-grade EEG device. Both signals are recorded simultaneously during 30 minutes long sessions in resting state. The data are collected from 22 patients with epileptiform abnormalities in EEG. In order to compare two multichannel EEG signals with differently placed references and electrodes, a novel data processing pipeline is proposed. It allows deriving matching pairs of time series which are suitable for similarity assessment through Pearson correlation. The average correlation of 0.55 is achieved on a test dataset, which can be considered a promising result, taking the positions shift due to the simultaneous electrode placement into account.

Introduction

Fast development of mobile and wearable sensors introduces new opportunities for field-trials and long-term monitoring in many research and clinical areas. However, these novel sensors require a rigorous validation, particularly when being used in clinical applications. The methods of validation can differ greatly depending on the type of the sensor and the purpose of the collected data. In this manuscript we focus on the validation of novel, wearable electroencephalography (EEG) sensors in a clinical environment.

In principle, two main approaches to **quantitative** EEG validation can be found in the literature:

Technical validation includes an assessment of electrical characteristics of the device or just the electrodes. This approach is often chosen by sensor developers or researchers as an initial validation. For example, Liao et al. (1) report the comparison of dry electrodes impedances to impedances of standard clinical wet electrodes. Generated signals of fixed frequency (2) can be used for a spectral analysis of the record. Play-back (1,3) is a technique where brain signal is recorded with the gold standard device and then re-played and recorded by the device in question. *In vivo* experiments are challenging, as the signal cannot be recorded at the same time and place with multiple devices. There are two possible solutions (3): “same time different place”, where several sensors are placed as close to each other as possible and data collected

simultaneously; and “same place different time”, where the recordings are being performed subsequently at the exact same positions. These methods are also referred to as parallel and serial (4). In the first case, the recorded signals can be compared through visual analysis (4,5), cross correlation (1,2,4–6), spectral features (7), mutual information (8,9) or other numerical features. In case of serial recordings, the options for quality assessment are more limited due to time-dependent changes in brain activity.

Another possibility is a **validation in a given experimental context**. It can be done to answer the question if the new device is suitable for a certain type of research or medical purposes. Here, the choice of possible tests is highly dependent on the application area, but typically the data are collected under well controlled conditions. Both parallel and serial setups can be used. For instance, the Emotiv Epoc mobile EEG has been investigated by Badcock et al. (10). The authors compare statistical features of auditory Event Related Potentials (ERP) captured in parallel experiments by the mobile and a research-grade device. Melnik et al. (2017) (11) used six different ERP paradigms and serial approach to study the variances caused by multiple sessions, subjects and devices. Another common option is an experiment when the subjects are instructed to have their eyes open and closed for a certain time interval. This allows comparing the EEG alpha band powers (2).

While all the above approaches are necessary and valid in a research and industrial environment, further validation is required prior to clinical use. Clinical validation proves that a device is capable of performing similar or better than an existing clinical gold standard under real-world conditions. Unfortunately, all described methods require performing the data collection under carefully controlled conditions, which is not always possible in a clinical environment. The present work aims to fill the methodological gap and provide a data analysis pipeline designed to compare the simultaneously recorded data from two multichannel EEG devices under the following limitations:

- No hardware modification can be done. In particular, the referencing electrodes are placed according to the respective design of the devices.
- No special external stimuli can be introduced, which limits statistical analysis of the signal responses.

The developed pipeline is investigated using simultaneously acquired data from a mobile and clinical EEG during a clinical study at the RWTH Aachen University Hospital.

Methods

Study design

Resting state EEG signals were collected simultaneously with a clinical and mobile EEG devices in a clinical study (12) from 22 patients, 11 females and 11 males. The patients were admitted to the hospital for epilepsy diagnostics. Average age of the patients was 40.2 ± 15 years. All patients were fluent in German language and signed the informed consent. The study was approved by the Ethical Board of Uniklinik RWTH Aachen (EK150-17 3.7.2017) and registered prospectively (DRKS-ID: DRKS00012424).

The clinical grade data were collected with Brain Quick Plus Evolution by Micromed, which is currently used as a standard device at Epileptology Section at Uniklinik RWTH (Aachen, Germany) where the recordings were conducted. This model has 21 EEG electrodes and reference electrode at G1. The electrodes are placed according to the 10-20 system (13) (see also Fig. 1a). The data are sampled at 256 Hz. 0.18 Hz high-pass filtering is performed in the amplifier. During patients' hospital stay they were undergoing continuous EEG monitoring for several days, with simultaneous collection of video data. Therefore we use the abbreviation vEEG for the clinical EEG further on. The video data were not used in our study. Approximately every 3 hours of EEG records are saved in a separate file and exported in the European Data Format (.EDF files).

For the mobile data collection, the Epoc Emotiv device was used (mEEG). This device was originally marketed as a mobile brain-computer interface (BCI) sensor, but gained a lot of interest from the research community (10,11). Despite a somewhat limited coverage of the cerebral cortex, the price range and number of electrodes (14, CMS/DRL references) made it a potentially useful tool for clinicians and researchers. The electrode positions are fixed due to a rigid plastic frame and approximate the 10-20 system (see Fig. 1a). Sampling rate of the device used in the study is 128 Hz.

In the experiment Epoc Emotiv was mounted while the clinical EEG device was already in use. mEEG was recorded for approximately 30 minutes within a single recording session. During these 30 minutes the patients were asked to limit their movements, as Emotiv Epoc device is highly prone to movement artifacts. At the beginning of the mobile data collection, the participants were asked to blink strongly several times in a row. This was done to facilitate the later data alignment.

The electrodes were placed as close as possible to the vEEG electrodes, at the same time minding the distance enforced by the glue used to attach vEEG electrodes.

The stiff frame of mEEG presented an additional challenge. After the mounting of mEEG a sketch of the placement was made (see Fig. 1a for an example).

Before any further processing, the data were pseudonymized and only the session's number was used to keep track of the files.

Three patients' data were discarded from the trial. One set because of the mobile recording software failure, one because of the failure of reference electrode in the clinical EEG, and the last one because of cell phone usage by the patient resulting in strong artifacts.

Data preprocessing

The raw EEG data files have different length, sampling rates and electrode placements and, therefore, require several preparation steps before the comparison between mobile and clinical signals can be done. The pipeline presented below allows for constructing aligned pairs of signal vectors, ready for further statistical analysis.

Initial pruning of clinical EEG data

The duration of the vEEG and mEEG signal is 3 hours and 30 minutes respectively. Based on the blinking and the marked time of the mEEG data collection, approximately 30 minutes of vEEG is cut (see Fig. 2). It is not yet perfectly aligned with the mEEG data, but the shift is within 15 sec.

Find corresponding electrodes between vEEG and mEEG.

The relative electrode positions of vEEG and mEEG are considered for each patient individually (see Fig. 1a). The couples of mobile and video electrodes placed directly next to each other (as close as the glue circle allows) are listed. For example, the following couples are chosen (Fig. 1b): mF3-vF3 and mT8-vT4. The set of two corresponding pairs of electrodes will be called a **quadruple**. The list of all possible electrode quadruples (e.g. [mF3, mT8, vF3, vT4], with m* and v* being electrodes from the mEEG and vEEG respectively) is created individually for each patient (Fig. 2).

Data cleaning and extraction of quadruples

Due to the relatively low signal-to-noise ratio, it is essential to reduce the noise before any comparison is made. Additionally, the data need to be synchronized in time and pruned to the same vector size. Therefore, the following data cleaning pipeline was implemented (see Fig. 3 for the visualisation):

- 1) The EEGLAB function *pop_rejchan* is used to detect corrupted channels in all files. The quadruples containing such channels are removed from the list.

- 2) Low pass filtering at $fq_lowpass$ Hz (the parameter $fq_lowpass$ depends on the chosen frequency band).
- 3) Downsampling vEEG to a sampling rate of 128 Hz, to match the lower mEEG sampling rate.
- 4) High pass filtering at $fq_highpass$ Hz (similarly to $fq_lowpass$, the parameter $fq_highpass$ depends on the chosen frequency band).
- 5) For each patient, the $perc$ -th percentiles of absolute values of the amplitudes are computed for mEEG and vEEG separately, and the corresponding raw data are divided by the resulting value. This normalization allows to avoid the problem of different scales of mEEG and vEEG.

As an output of the described process, we obtain filtered, normalized and roughly time-aligned mobile and clinical EEG data sets for each patient. Additionally the corrupted channels are eliminated and an individual list of channel quadruples to be used for further re-referencing is stored.

Processing of single quadruples

More precise time alignment and artifact removal is done for individual quadruples (e.g. [mF3, mT8, vF3, vT4] as illustrated on Fig. 1b) of specific patients. The resulting vectors can be compared through Pearson correlation.

Artifact detection for a single channel

In order to capture short-term signal disturbances, we construct an adjusted procedure for artifact detection:

- 1) First, threshold **AmThresh** is fixed and all data with absolute value exceeding AmThresh is marked as NaN (not a number). The data points are not removed to allow later time synchronization. It should be noted, that previous data normalization makes it possible to choose for a common single threshold for both mEEG and vEEG records.
- 2) Next, non-overlapping intervals of length **WinLength** are taken, average of amplitudes' absolute values is computed and the whole window is marked as NaN if this average exceeds the parameter **AmThreshWin**.
- 3) **Artifact index** is computed by dividing the length of the "corrupted" data (marked as NaN by the total data vector length.
- 4) If **Artifact index** exceeds 70%, all quadruples including this channel are removed from the list.

The described procedure is subsequently applied for all four time series involved in the quadruple.

Remark 1. One of the most common procedures for EEG analysis is the removal of eye blinking artifacts. Here, it was decided against this removal, because the locations typically known for

strong eye artifacts usually are not involved in the analysis. Additionally, such artifacts are a normal part of the EEG signal, and as such should manifest similarly in mEEG and vEEG signals. One may argue, that high amplitude of such signals may unproportionally influence the linear correlation, but since they are not very prominent in this particular data they were neglected. Similarly, ECG artifacts were not significant in the considered data.

Re-reference of quadruples.

We mathematically re-reference the electrodes within each quadruple (through pairwise subtraction) to eliminate the effect of the global reference (e.g. quadruple [mF3, mT8, vF3, vT4] equivalent to [mF3-mRef, mT8-mRef, vF3-vRef, vT4-vRef] after re-referencing results in a pair [mF3-vF3, mT8-vT4]). This subtraction can be thought of as bipolar (BP) re-referencing. Thereby, the global references are eliminated. The spatial proximities within the quadruple should result in a similarity of such differences. Therefore, the time series mF3-mT8 is expected to show strong similarities with time series vF3-vT4. Note, that both time series may contain NaN terms propagated from the artifact removal procedure. From here on, mEEG and vEEG will still refer to mobile and clinical EEG data in the original referencing while mEEG(BP) and vEEG(BP) will refer to respective data re-referenced in bipolar (BP) manner. Note that the bipolar referencing is different for each patient and is based on patient-specific relative electrode placement. Similarly, the number of quadruples may vary.

Fine aligning

At this stage, MATLAB function *lag* can be used to find a time shift between the video and mobile data. This shift is assumed to remain constant throughout the session.

Remark 2. Finding the correct lag is only possible if the data has a certain minimal signal-to-noise ratio, otherwise all correlations are close to zero and the time shift is set to zero by the algorithm. For one recording session it is possible that the lags slightly vary between quadruples, but typically the differences are within 1-2 data points. Therefore the lags are first computed for all quadruples of a given patient, then the majority vote is used to set the same time shift to all combinations, allowing a fluctuation of ± 5 data points to compensate for differences in device-specific recording order of the channels.

Remark 3. If the frequency band chosen for filtering is too narrow, errors might occur in the lag computation. Therefore it is recommended to pre-compute and save lags for a wider band.

Finally, the longer vEEG vectors are pruned to match the length of the mEEG vectors and the time stamps from the mEEG data are added. NaN data points from each single time series in the quadruple propagated to the bipolar re-referencing. After the aligning, all data segments where either of the two bipolar time series contained NaN were removed (see Fig. 4).

Finally, Pearson linear correlation for the given quadruple is computed and stored.

Aggregating data for multiple patients

For each patient and fixed set of the parameters, the table containing all quadruples and their corresponding correlation coefficients is built. The patient average and the grand average across all patients can be reported. Additionally, the data loss resulting from all processing steps is tracked.

Optimization and parameters choice

Let us notice, that there are 6 parameters which we can choose in the above described procedure: *fq_lowpass*, *fq_highpass*, *perc*, *AmThresh*, *WinLength* and *AmThreshWin*. A procedure for establishing parameters values consists of several steps.

First, we fix the default filtering parameters *fq_lowpass* and *fq_highpass* based on the estimation of the usable spectrum of the data. In principle the frequency parameters can be chosen freely, depending on which part of the spectrum the correlation is of interest. For instance, we will argue later, that alpha band (7.5-12.5 Hz) shows the best correlation. Nevertheless, in order to find the broadest default frequency band reasonable for our research, the average cross-spectrum of the data is computed and thresholding is made to fix the initial *fq_lowpass* and *fq_highpass* values (Fig. 5). This cross-spectrum is computed in a procedure similar to the above described correlation computation, with only the initial high-pass filtering of mEEG at 0.5 Hz. Other filtering, normalization and artifact removal are skipped for the estimation of the filter values. Second, a set of discrete values for *perc*, *AmThresh*, *WinLength*, *AmThreshWin* is chosen. For each parameter combination and patient, the correlation tables are computed.

Third, the data are randomly divided into two subsets of 14 training and 5 test records. The parameter set with the best resulting grand average correlation coefficient is established based on the first set and quasi-independently evaluated based on the second set.

Correlations between vEEG electrodes (vEEG to vEEG correlation)

As direct performance comparison, we also investigate the similarities within the vEEG signal and calculate the correlations. In this case, all neighbouring pairs of vEEG electrodes were defined (e.g. [vP4, vT6]), and underwent regular artifact removal procedures. Instead of quadruples that are re-referenced between two electrodes from both vEEG and mEEG each, we use two vEEG electrodes and their reference electrode. Since both electrodes are referenced to the same reference electrode, this constructs a virtual quadruple (e.g., [G1-vP4, G1- vT6] with G1 being the reference electrode) and correlation between the two signals can be calculated.

In case of vEEG-mEEG, the difference between the computed correlation and perfect linear correlation value (1) is due to the following factors: a) devices quality differences, b) noise, c) spatial distances between two pairs within one quadruple (e.g. vF3 to mF3 and vT4 to mT8). In case of vEEG-vEEG correlation, the differences are due to a) noise and b) spatial distances

between the two vEEG electrodes (e.g. vP4 to vT6). These distances are larger than in vEEG to mEEG pairs, but there is no spatial distance influencing the correlation for the reference electrode. If we assume that these differences give a comparable error and the noises are on the same level, then the differences in average correlation should reflect the differences in the devices quality.

Software

Data collection was performed with BrainLab (14). Data analysis was done with MATLAB R2017b (The Mathworks, Natick, MA, USA) and EEGLAB v.14 (15).

Results

Establishing frequency interval for preprocessing

In the first step the frequency band with the boundaries $fq_highpass$ and $fq_lowpass$ for the preprocessing is chosen for the analysis. The lower boundary ($fq_highpass$) is set to 1Hz, which is sufficient to remove the drift from the EEG data. In order to set $fq_lowpass$, the magnitude-squared coherence was calculated to assess which frequencies are useful for comparison (Fig. 5). The magnitude-squared coherence indicates the shared information between mEEG and vEEG in distinct frequency bands. A cut-off value of 0.1 was chosen, resulting in upper frequency band limit of 38Hz. This is in line with the commonly used frequency bands in EEG analysis (following (21), delta up to 3.5Hz, theta 3.5-7.5Hz, alpha 7.5-12.5Hz, beta 12.5-30Hz and gamma > 30Hz). Recent research also shows the importance of higher gamma (> 50Hz) in certain application contexts (16,17), but the chosen interval of 1-38Hz is sufficient for the basic neurologic assessment.

Parameter optimization

For the chosen frequency band (1-38Hz), we consider the following discrete values of the parameters:

$perc = \{80\%, 85\%, 90\%, 95\%, 97\%, 99\%\}$,
 $AmThresh = \{1, 1.5, 2.0, 3.0, 4.0, 5.0, 6.0, 7.0\}$,
 $AmThreshWin = \{0.5, 1.0, 1.2, 1.5, 2, 2.5\}$,
 $WinLength = \{15, 50, 75, 100, 200, 250\}$.

For all possible parameter combinations, the average correlation was computed on the randomly chosen train set. The following optimal set of parameters, resulting in mean correlation of 0.46 was obtained through exhaustive search:

perc= 85%,
AmThresh=6.0,
WinLength=100,
AmThreshWin=1.5.

As a boundary condition, the amount of data retained needed to be 75% on the train data.

Correlation between vEEG and mEEG

The optimal parameters determined on the train set were fixed for evaluation on the test set. On the test set, the resulting mean correlation equals to 0.55. The volume of data left for the analysis after artifact removal equals 76% of the original data volume.

In order to better understand the distribution of correlations, each patient was analyzed individually (Fig. 6). While most patients produce similar results, some patients have much lower overall correlation coefficients.

The results are highly sensitive to the choice of the frequency band. Considering the fact, that the bands of interest might differ depending on the application, the correlations were examined independently for each frequency band (Tab. 1). The alpha band shows the highest correlation.

Correlations between vEEG electrodes

In order to obtain the comparison correlation values, the correlation between neighboring video electrodes were computed and averaged. The received grand average was 0.57 (see Fig. 7 for details). To compare the quality of the mEEG and vEEG, we report the percentage of data retained during cleaning for both mEEG to vEEG and vEEG to vEEG correlations (Tab. 2). The vEEG to vEEG comparison retains approx. 88-95% of the data, which means that 5-12% are discarded. In the mEEG to vEEG comparison, approx. 25% of the data are discarded, which indicates that the mEEG is responsible for about 15% or two thirds of the total discarded data.

Discussion

Main Findings

The main goal of this work was to find a way to objectively (quantitatively) validate a mobile EEG device in a situation when only resting state brain waves can be collected and no hardware manipulations can be performed (e.g., due to regulations). In particular, the positions of the

individual electrodes are fixed and no common referencing is available. To overcome these limitations, we have developed an advanced and multi-parametric data processing procedure, which allows obtaining illustrative and robust results in the form of a grand average correlation coefficient.

The resulting average correlation is dependent on the chosen frequency range. In case of EEG data it is reasonable to consider standard brain wave frequency bands (delta, theta, alpha and beta). The average correlation varies from 0.52 on the delta band of the test set to 0.66 on alpha band. The average across the overall considered frequency band (1-38Hz) is 0.55. While this number may look moderate, it can be perceived as quite high if we take into account that as much as 76% of all data are preserved and that we compute a grand average over multiple (on average 53 per patient) pairs of 30 min long data vectors (test set). Due to chance, the average correlation on the test set is higher than on the train set. One of the main reasons is a consistency in signal quality within one record and relatively small ($n=19$) number of records, resulting in only 5 test records. Nevertheless, both train and test sets show similarity in the results and have similar dynamics across the bands, with alpha band showing the highest correlation.

Although the overall data quality was not optimal, only three out of 22 sessions were rejected because of recording failures, but none during the processing pipeline execution. For the remaining 19 sessions, on average, 53 quadruples of electrodes were selected based on strictly formulated criteria. Furthermore, only 24-28% of data were lost due to rejection of full channels or signal segments (train and test sets). The rejection was performed only based on computer algorithms. Thereby, the reported correlation averages across multiple session and quadruples, and covers most of the recorded data. With respect to the high amount of data and potential variation, the low standard deviation of 0.19 (on test set) points towards high robustness.

Due to the uniqueness of the data it is a difficult task to evaluate the results. Therefore we have performed a test to relate the grand average correlation. The average correlation was computed for neighbouring pairs of vEEG electrodes resulted in a correlation value of 0.57, which is comparable to vEEG-mEEG correlation. The assumptions regarding possible deviations from a perfect correlation we have listed above are difficult to verify, so we can only hypothesize that these close values suggest comparable signal quality of mobile and clinical devices.

Comparison to the state of the art

Analyzing other works on electrodes quality comparison, where experiments were performed in a more controlled environment, may bring a better understanding of how much of the signal differences are caused by spatial shift of the electrodes and how much by the lower quality of the mobile hardware. Only experiments of “same-time-different-place” type with resting state or similar conditions were chosen for the comparison. In some papers, the authors compare the same type of electrodes in order to estimate the spatial shift-related signal change.

In Fiedler et al. (6) the authors place 3 types of dry electrodes at Fp1, Fp2, O1 and O2 sites, with standard wet electrodes adjacent. Additionally two sets of wet electrodes was tested to provide a baseline. In resting state the resulting average correlations were: 0.24, 0.59, and 0.25 for three dry-to-wet comparisons respectively and 0.58 for wet-to-wet combination.

Estepp et al. (5) under the open eyes condition reported 0.84, 0.61 and 0.32 for dry-to-wet comparison at Fz, C4 and Pz positions respectively. Similarly, wet-to-wet combinations resulted in 0.97, 0.95 and 0.80.

In Wyckoff et al. (2) the measurements with dry and wet electrodes were done at Fz, C3, Cz, C4 and Pz. For the open eyes condition the average correlation varied from 0.28 on delta band to 0.99 on alpha, beta 1 (13-16Hz) and beta 2 (13-21Hz) bands.

Liao et al. (1) reported respectively 0.95 and 0.91 correlation at F10 and POz for two different electrode types.

More references can be found in a review paper of Lopez-Gordon et al. (3). The results in the above referenced papers are characterized by high variability. The reported correlations vary from 0.25 to 0.97 on time-domain signal, which can be explained by different quality of the tested electrodes, but also by differences in the placement, experimental details and data processing. All the experiments were carefully controlled, electrode number limited to a maximum of 5, subjects movement could be minimized and sometimes the segments of data rejected after visual examination (4).

In contrast, in our research, the data were collected under minimal control, multiple electrode sites were used (including the ones known for high artifact presence) and no human examination was used for data processing. Yet, the results are revealing a correlation level comparable to other comparisons with a similar setup.

Limitations

In the described work Pearson correlation was used as a straightforward similarity measurement. It seems to be a reasonable choice in the given context, as linear relationship is exactly what we are expecting from our experiments if the spatial shift is neglected. On the other hand, this spatial shift may introduce significant nonlinear effects. In the further work it might be beneficial to consider different measurements of similarity, such as mutual information (9). Nonlinear relationships in EEG have been deeply studied in the context of epileptic seizures where synchronisation of different brain areas often occurs. In Quiroga et al. (8) a number of nonlinear measurements is discussed, however the same paper also suggests ultimate resulting similarity of the different types of measurements, including linear correlation.

Another limitation is the restriction of the comparison bandwidth from 1 to 38Hz. The maximum available frequency of our setup is 64, however only little mutual information was detectable in the higher frequencies. This indicates both a limitation in the hardware, as well as in the proposed algorithm, as the parts with potentially worst correlation are excluded. While this is reasonable as long as the values chosen are still clinically relevant, it requires special care when choosing the *fq_lowpass* value.

Conclusion

In this work we have developed a data analysis procedure designed to deal with two sources of EEG data recorded simultaneously during resting state. This procedure aims to provide an objective measurement of the data quality.

It is not uncommon that in the procedure of EEG comparison visual assessment by the trained specialists is used as a part of data pre-processing (4). While professional opinion may provide a unique insight, it is also very costly to obtain. In the particular case of our study, more than 20 hours of multichannel data need to be analyzed. Moreover, multiple studies have shown, that human assessment is not fully reproducible and high intra- and inter-rater variances are consistently reported (18,19), which could be reproduced on our data in a previous publication (12). In contrast, automated analysis pipeline provides an objective, fast and low-cost way to perform the data comparison.

The presented procedure deals with the challenges of different referencing, spatial shifting of the electrodes and lack of controlled stimuli (such as in ERP experiments). Using automatically optimized parameters for the pre-processing, a grand average linear correlation of 0.55 between mobile and clinical EEG devices was obtained. It was compared to several baseline correlations, such as clinical-to-clinical EEG correlation, to conclude that the overall quality of the considered mobile device is good, since similar correlations can be seen when only electrode types are changed. This result agrees with our previous study, where trained neurologists were clinically investigating the data (12), and with a number of other studies done on Epoc Emotiv in different contexts (11,20). However, to our best knowledge, a fully automated approach in combination with the resting state data was not previously reported.

The presented pipeline might benefit in the future from including more sophisticated signal processing methods, such as mutual information. Nevertheless, in the current form it already shows high efficiency and can be easily generalized to different multi-channel sensors, such as EMG or ECG.

Acknowledgements

The authors are grateful to Jonas Schulte-Coerne for valuable discussions on the signal processing.

References

1. Liao L-D, Wang I-J, Chen S-F, Chang J-Y, Lin C-T. Design, Fabrication and Experimental Validation of a Novel Dry-Contact Sensor for Measuring Electroencephalography Signals without Skin Preparation. *Sensors*. 2011 May 30;11(6):5819–34.
2. Wyckoff SN, Sherlin LH, Ford NL, Dalke D. Validation of a wireless dry electrode system for electroencephalography. *J NeuroEngineering Rehabil* [Internet]. 2015 Oct 31 [cited 2018 May 25];12. Available from: <https://www.ncbi.nlm.nih.gov/pmc/articles/PMC4628242/>
3. Lopez-Gordo M, Sanchez-Morillo D, Valle F. Dry EEG Electrodes. *Sensors*. 2014 Jul 18;14(7):12847–70.
4. Gargiulo G, Calvo RA, Bifulco P, Cesarelli M, Jin C, Mohamed A, et al. A new EEG recording system for passive dry electrodes. *Clin Neurophysiol*. 2010 May;121(5):686–93.
5. Estep JR, Christensen JC, Monnin JW, Davis IM, Wilson GF. Validation of a Dry Electrode System for EEG. *Hum Factors Ergon Soc Annu Meet Proc*. 2009 Oct 19;53(18):1171–5.
6. Fiedler P, Hauelsen J, Jannek D, Griebel S, Zentner L, Vaz F, et al. Comparison of three types of dry electrodes for electroencephalography. *ACTA IMEKO*. 2014 Sep 23;3(3):33.
7. Grozea C, Voinescu CD, Fazli S. Bristle-sensors—low-cost flexible passive dry EEG electrodes for neurofeedback and BCI applications. *J Neural Eng*. 2011;8(2):025008.
8. Quian Quiroga R, Kraskov A, Kreuz T, Grassberger P. Performance of different synchronization measures in real data: a case study on electroencephalographic signals. *Phys Rev E Stat Nonlin Soft Matter Phys*. 2002 Apr;65(4 Pt 1):041903.
9. Mikkelsen KB, Kidmose P, Hansen LK. On the Keyhole Hypothesis: High Mutual Information between Ear and Scalp EEG. *Front Hum Neurosci* [Internet]. 2017 Jun 30 [cited 2018 May 25];11. Available from: <http://journal.frontiersin.org/article/10.3389/fnhum.2017.00341/full>
10. Badcock NA, Preece KA, de Wit B, Glenn K, Fieder N, Thie J, et al. Validation of the Emotiv EPOC EEG system for research quality auditory event-related potentials in children. *PeerJ*. 2015;3:e907.
11. Melnik A, Legkov P, Izdebski K, Kärcher SM, Hairston WD, Ferris DP, et al. Systems, Subjects, Sessions: To What Extent Do These Factors Influence EEG Data? *Front Hum Neurosci* [Internet]. 2017 Mar 30;11. Available from: <http://www.ncbi.nlm.nih.gov/pmc/articles/PMC5371608/>
12. Titgemeyer Y. Can commercially available wearable EEG devices be used for diagnostic purposes? An explorative pilot study.
13. Tatum WO. *Handbook of EEG interpretation* [Internet]. New York [N.Y.: Demos Medical Pub.; 2013 [cited 2018 Oct 8]. Available from: <http://www.credoreference.com/book/spheegi>
14. Fink I. *Mobile EEG Presentation Framework* [Bachelor Thesis]. [Aachen, Germany]: RWTH Aachen University; 2016.
15. Delorme A, Makeig S. EEGLAB: an open source toolbox for analysis of single-trial EEG dynamics including independent component analysis. *J Neurosci Methods*. 2004 Mar 15;134(1):9–21.
16. Ball T, Demandt E, Mutschler I, Neitzel E, Mehring C, Vogt K, et al. Movement related activity in the high gamma range of the human EEG. *NeuroImage*. 2008 Jun;41(2):302–10.
17. Darvas F, Scherer R, Ojemann JG, Rao RP, Miller KJ, Sorensen LB. High gamma mapping using EEG. *NeuroImage*. 2010 Jan 1;49(1):930–8.

18. Grant AC, Abdel-Baki SG, Weedon J, Arnedo V, Chari G, Koziorynska E, et al. EEG interpretation reliability and interpreter confidence: a large single-center study. *Epilepsy Behav* EB. 2014 Mar;32:102–7.
19. Benbadis SR. Generalized Epilepsies on EEG: Overview, Types of Waveforms Associated With Generalized Epilepsies, Clinical Presentation. 2017 Mar 30 [cited 2017 Apr 24]; Available from: <http://emedicine.medscape.com/article/1140724-overview>
20. McKenzie ED, Lim ASP, Leung ECW, Cole AJ, Lam AD, Eloyan A, et al. Validation of a smartphone-based EEG among people with epilepsy: A prospective study. *Sci Rep*. 2017 Apr 3;7:srep45567.
21. Zschocke S, Kursawe H, editors. *Klinische Elektroenzephalographie*, Springer Medizin; 2012. 656 p.

Figure 1

Relative electrode positions.

(A) Example of the vEEG (green) and mEEG (blue) electrode placement. The placements differ slightly from patient to patient, therefore a sketch is made for each session to track the relative electrode positioning. (B) Zoom to an example electrode quadruple.

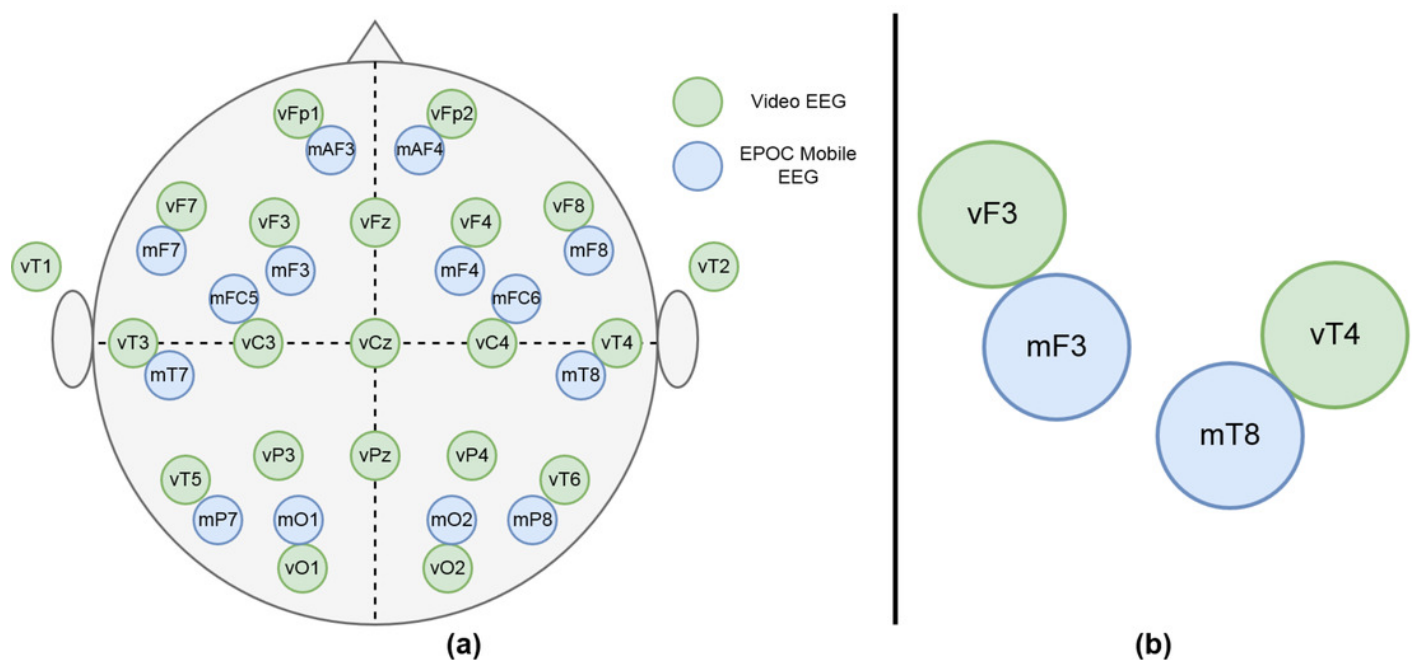


Figure 2

Initial pruning and quadruplet list generation.

Using timestamps and blinking artifacts approximately simultaneous vEEG and mEEG files are obtained. Personal sketch of electrode positions allows to detect spatially close electrode pairs and generate reference-free quadruples.

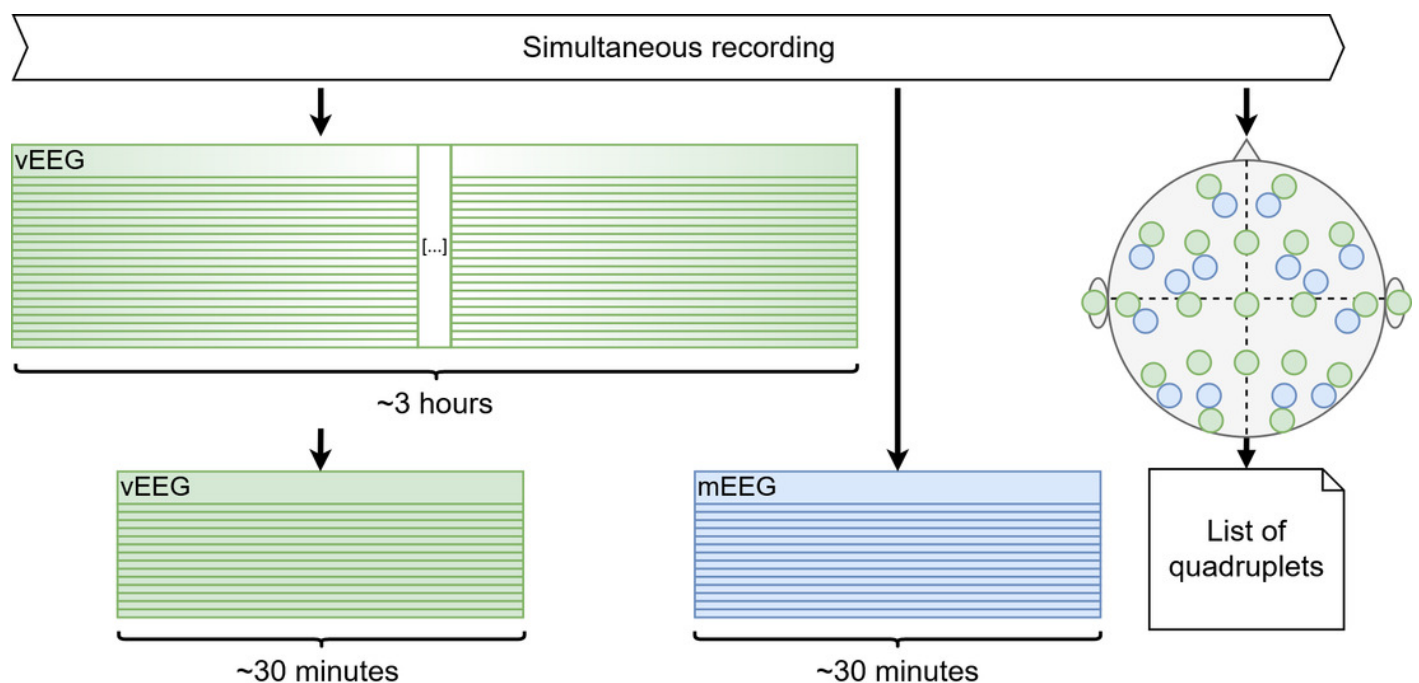


Figure 3

Data cleaning.

Data cleaning involves removal of corrupted channels (followed by the update of the quadruple list), filtering, re-sampling and normalization of the data.

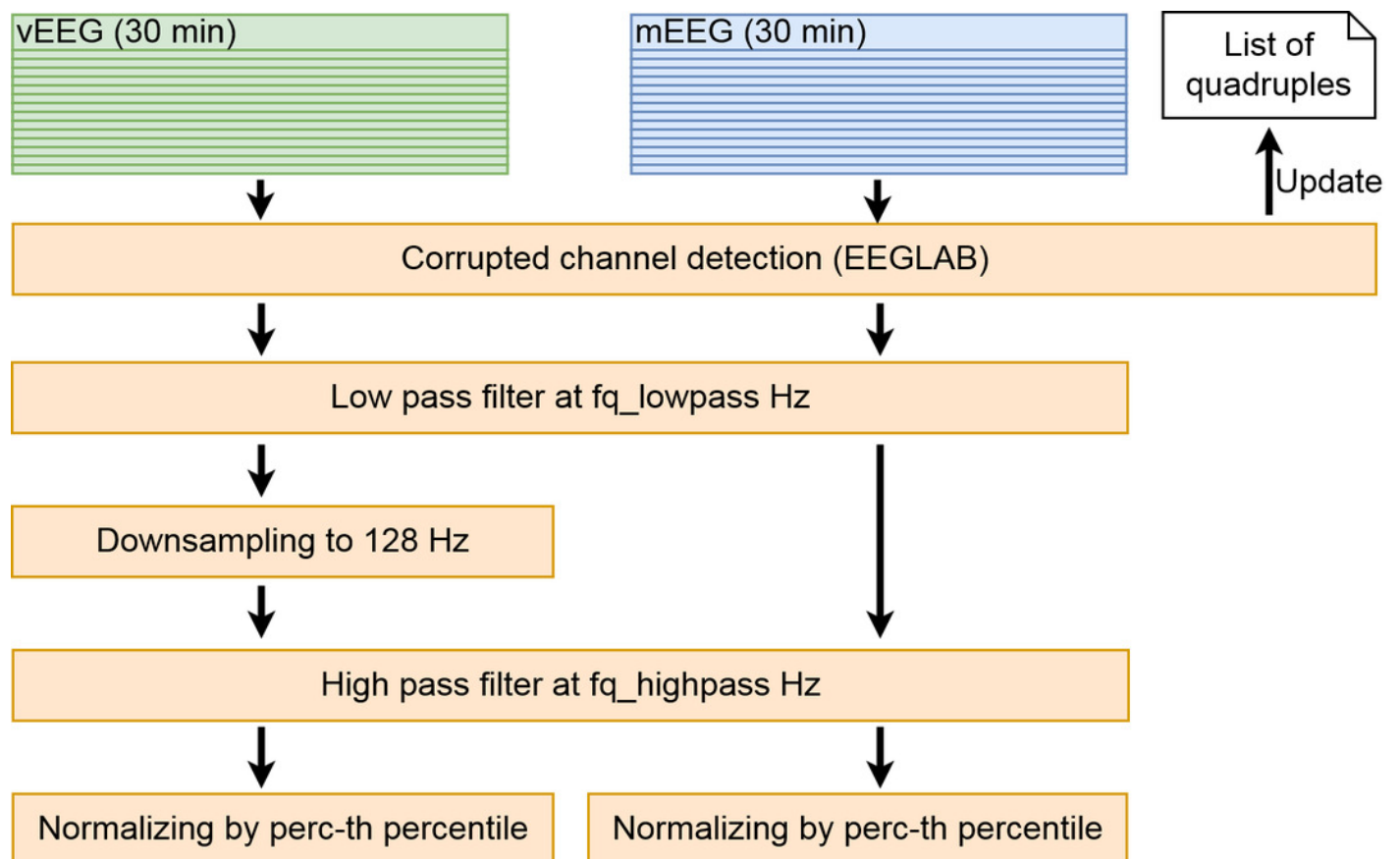


Figure 4

Pipeline for creating matching data vectors to compare through linear correlation.

For each quadruple artifact detection is performed in each electrode's data, then time aligned reference free vectors are constructed and the artifacts are removed without breaking the time synchronisation. Vectors with low signal-to-noise ratio are removed from the list. Linear correlation is computed for the remaining matching vectors.

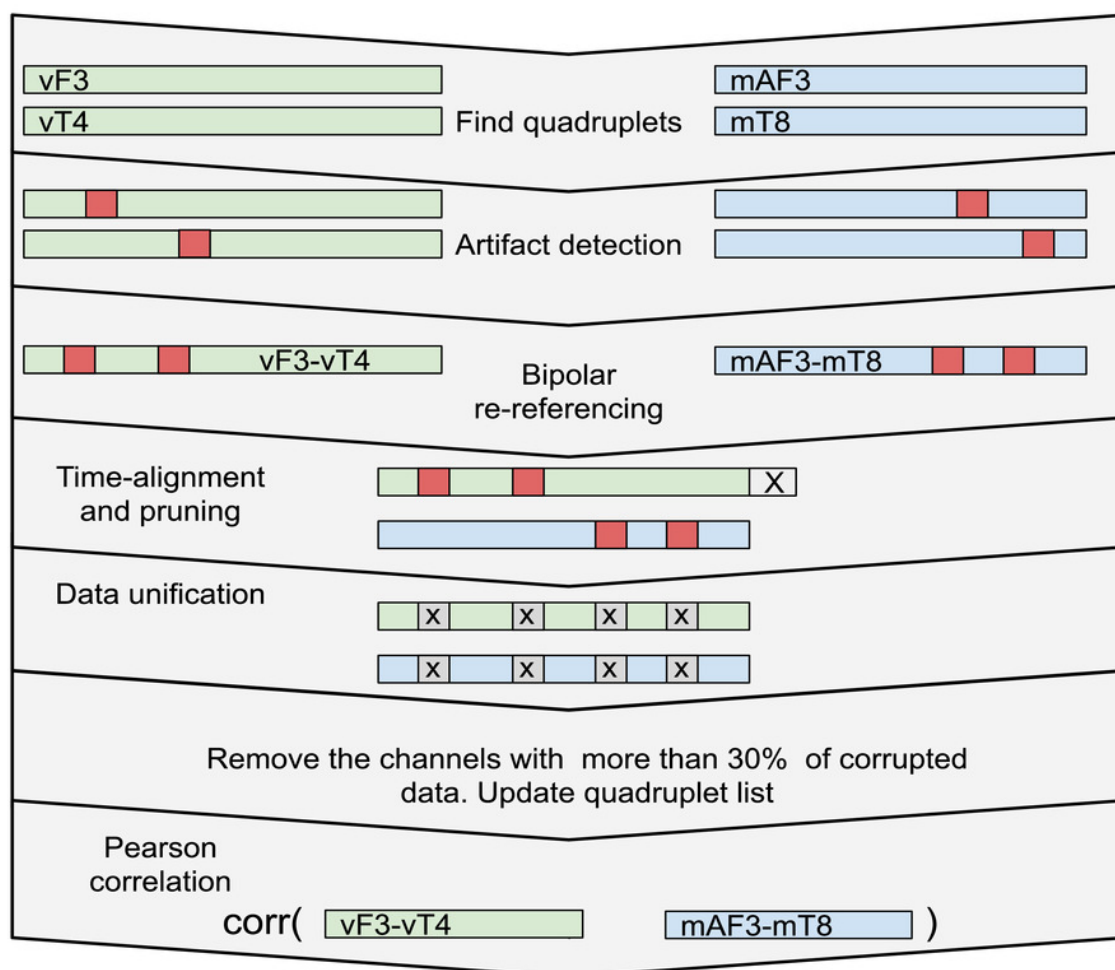


Figure 5

Magnitude-squared coherence between mEEG and vEEG averaged across electrodes and patients.

Red lines indicate thresholds for $f_{q_highpass}$ and $f_{q_lowpass}$.

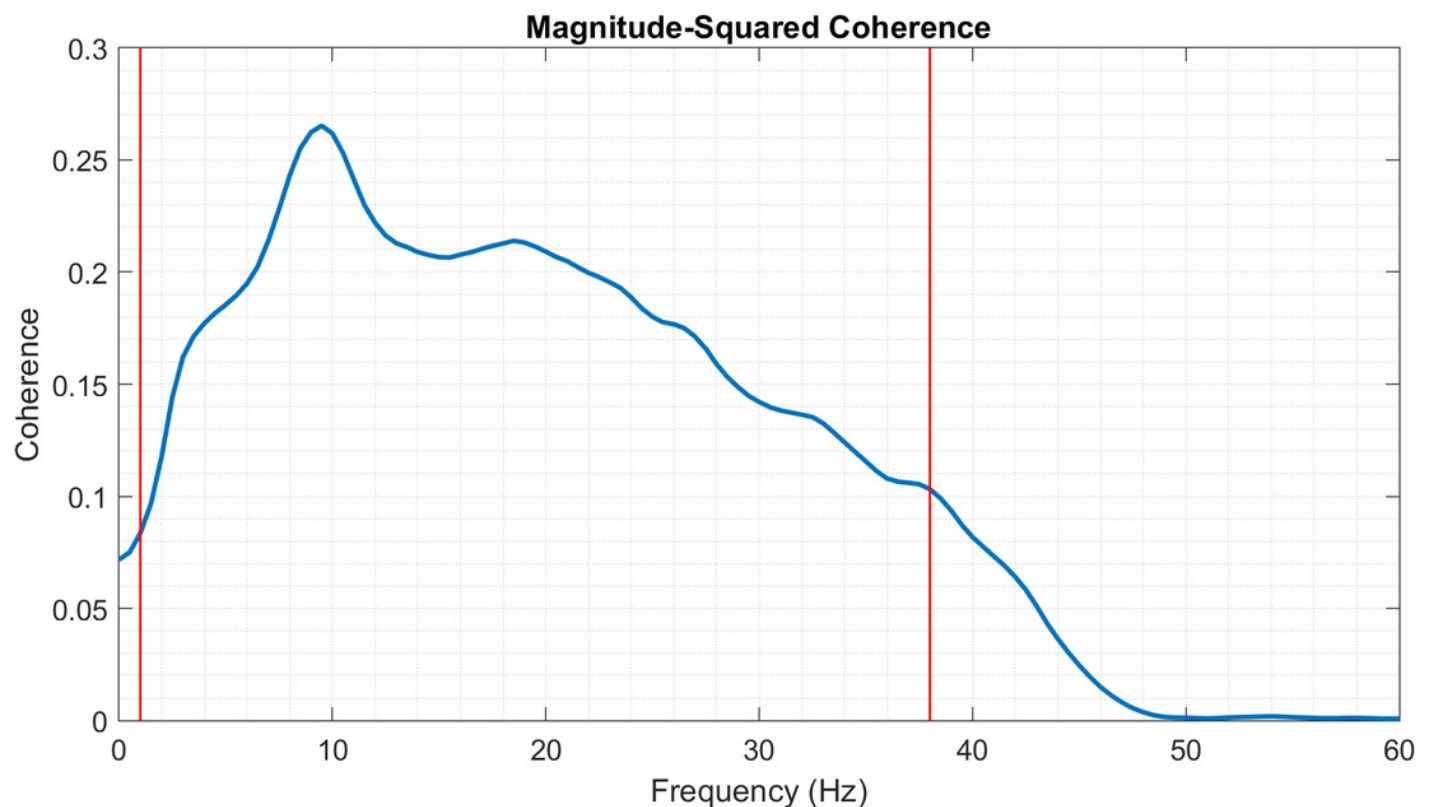


Figure 6

Distribution of correlation for each quadruple per individual patient.

Results on the test set are shown in red and on the train set in blue.

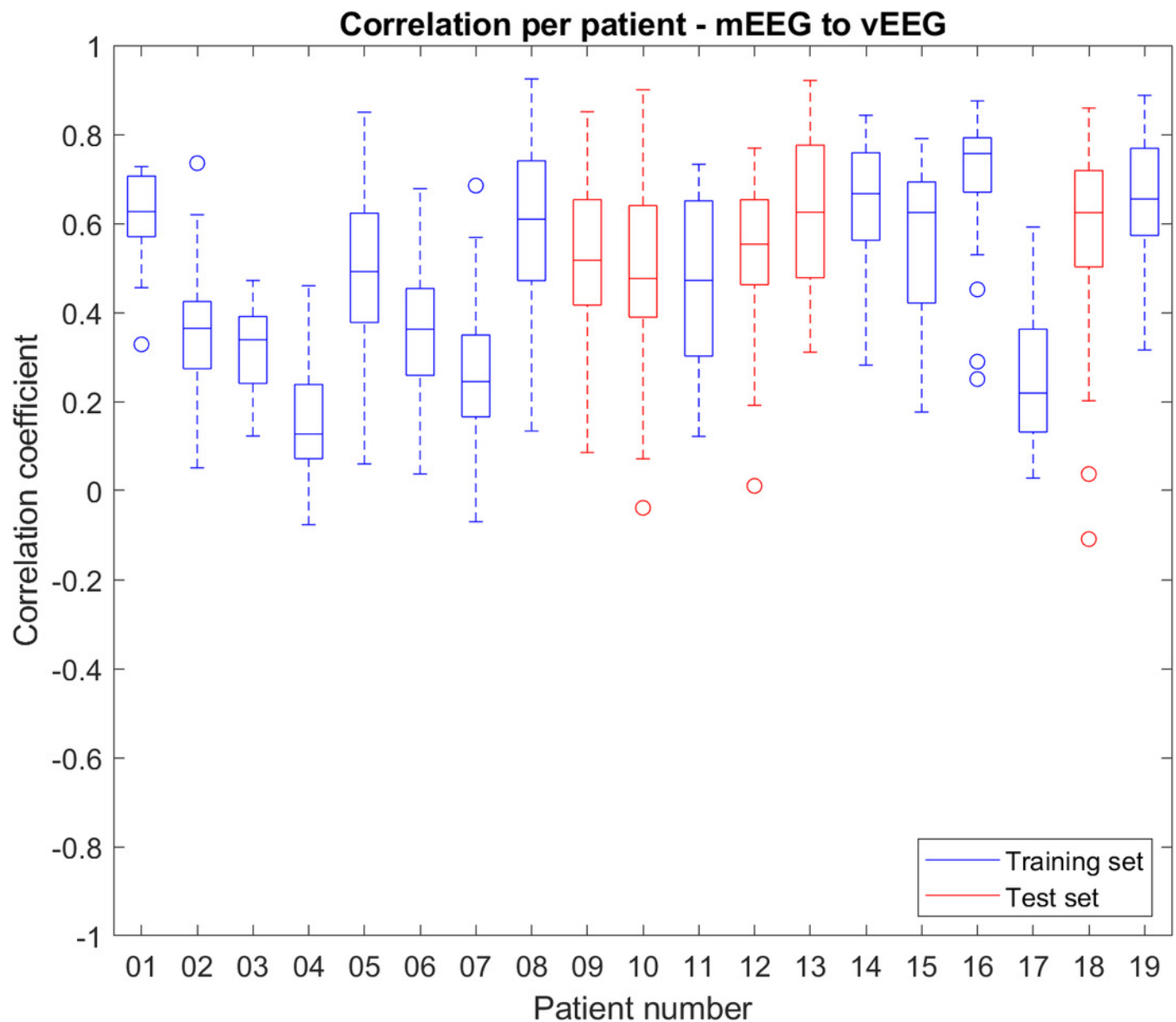


Figure 7

Distribution of correlations over individual patients.

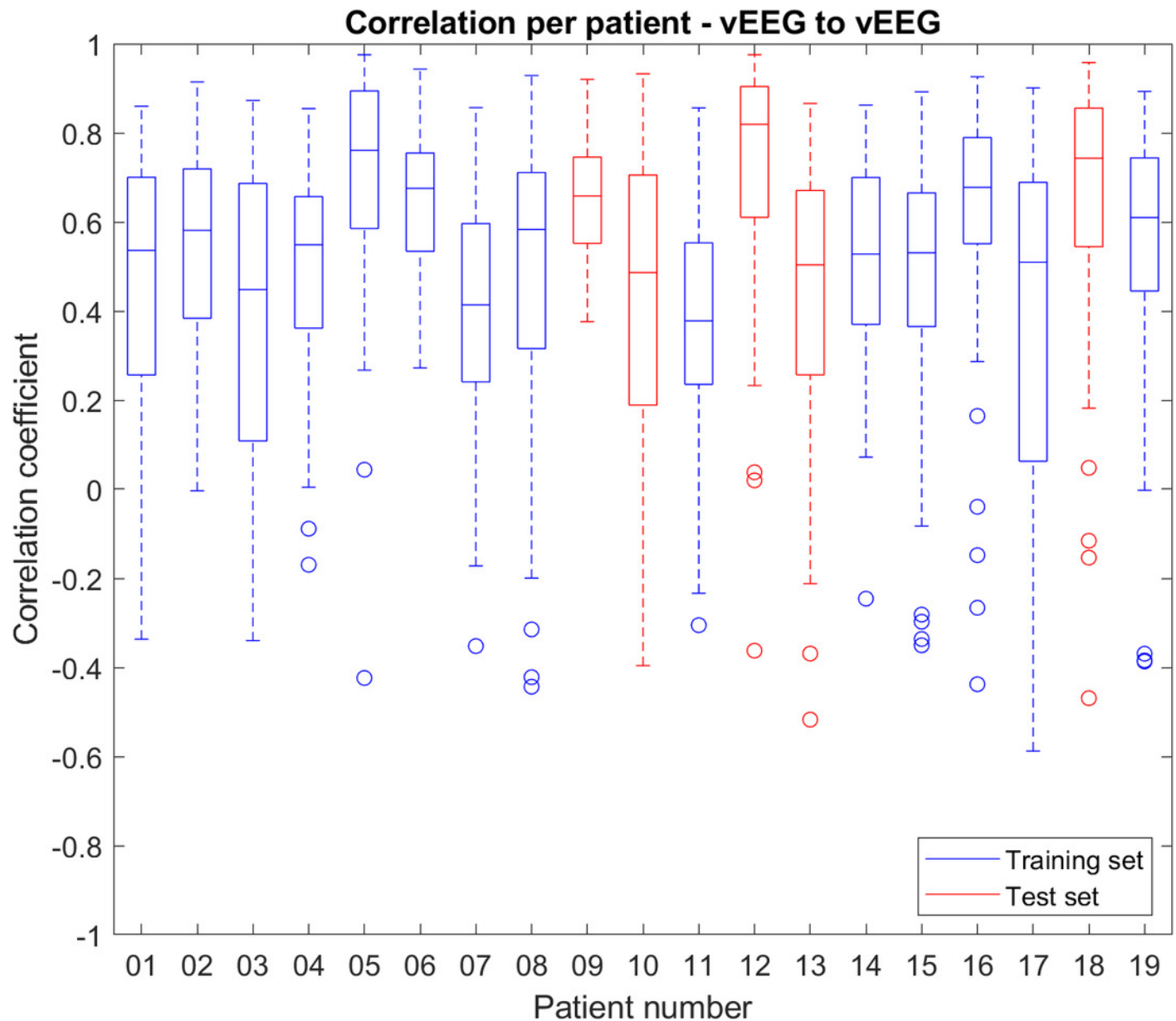


Table 1(on next page)

Average correlation by frequency band for the test set group. Gamma is not analyzed individually as only part of the band is present in the signal.

Average correlation	Delta (1-3.5Hz)	Theta (3.5-7.5Hz)	Alpha (7.5-12.5Hz)	Beta (12.5-30Hz)	Full band (1-38Hz)
Train set	0.40	0.51	0.56	0.45	0.46
Test set	0.52	0.63	0.66	0.57	0.55
Full data set	0.43	0.54	0.59	0.49	0.49

1
2

Table 2(on next page)

Percentage of data left after channel rejection and artifact removal.

Data left in total (in %)	Test set	Train set	Full data set
mEEG - vEEG	75.80	72.17	73.13
vEEG - vEEG	95.25	88.65	90.38

1
2